



GASTROSTOMY CATHETER

100% Silicone for Enhanced Patient Comfort



GASTROENTEROLOGIC
SERIES

GASTROSTOMY CATHETER

Fortune® Silicone Gastrostomy Catheter:


- 100% medical grade silicone for superior biocompatibility.
- Proper softness and excellent flexibility for enhanced patient comfort.
- Transparent tube allows for easy visual inspection. Radiopaque allows for confirmation the position of tube.
- Depth marking for identification of the tube position.
- Wide open-rounded end improves flow and minimizes trauma on insertion.
- Shortened tip design reduces gastric mucosa irritation.
- Funnel with several inner sealing o-ring for other device connecting universally.
- Attached plug for secured sealing to avoid air entering and contamination.
- External Bolster prevents tube migration and eliminates the needs for tape or suture.
- It is not easily visible when wearing normal clothes and is more comfortable.
- Universal medication port.
- Latex/DEHP free.



SPECIFICATIONS

Gastrostomy Catheter:

REF. No.	Size	Balloon volume	Length	Depth marking	X-RAY	Bolster
2019-0012	12FR	5-10cc				
2019-0014	14FR					
2019-0016	16FR	10-20cc		2		
2019-0018	18FR			4		
2019-0020	20FR		255 mm	6	yes	yes
2019-0022	22FR	15-20cc		8		
2019-0024	24FR			10		
2019-0026	26FR			cm		
2019-0028	28FR					

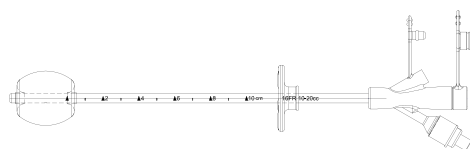

2019-00XX



6 F., No. 29, Sec. 2, Jhongjheng E. Rd.,
 Danshuei Dist, New Taipei City, 251 Taiwan
 Tel.:(886) 2-2624-2233 Fax.:(886) 2-2624-2266
 www.fortunemed.com Mail:info@fortunemed.com

GASTROSTOMY CATHETER WITH ENFIT CONNECTOR

Fortune® Silicone Gastrostomy Catheter with ENFit

Connector:


- The silicone gastrostomy catheter with purple ENFit connector is soft and flexible, It includes a silicone bolster, a silicone balloon, and silicone tri-ports connector. These components are medical grade silicone for enhancing patient comfort and superior biocompatibility; It can reduce the discomfort contact with the tube.
- Transparent tube for easy visual inspection, X-ray-opaque line for easy and accurate confirmation of tube position by X-ray. Depth markings for convenient confirmation of the tube position.
- Wide and open-rounded end improves flow, and the shortened tip design reduces gastric mucosa irritation.
- The tri-ports connector of the gastrostomy tube has assembly 2 ENFit male connectors and 2 ENFit female caps in the feeding lumen and the medication lumen. When the feeding or medication is stopped, the ENFit female caps can be close by rotated clockwise to prevent the air into the stomach.
- External Bolster prevents tube migration and eliminates the needs for tape or suture.
- The ENFit connector can prevent the misconnection risk, and provide a safer feed and medication channel for the patient.



SPECIFICATIONS

Gastrostomy Catheter with ENFit Connector:

REF. No.	Size	Balloon volume	Length	Depth marking	Description
2019-0512	12FR	5-10cc			
2019-0514	14FR				- Silicone Bolster
2019-0516	16FR	10-20cc		2	- ENFit Male Connector*2
2019-0518	18FR			4	- ENFit Female Cap*2
2019-0520	20FR		256 mm	6	- FEED Feeding
2019-0522	22FR	15-20cc		8	- MED Medication
2019-0524	24FR			10	- Depth Marking
2019-0526	26FR			cm	- X-Ray Opaque Line
2019-0528	28FR				

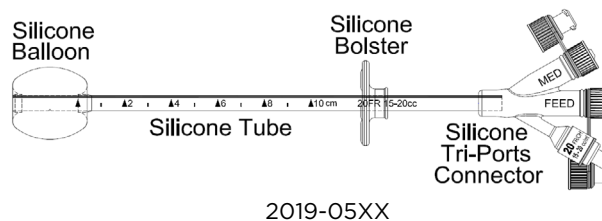

6 F., No. 29, Sec. 2, Jhongjheng E. Rd.,
Danshuei Dist, New Taipei City. 251 Taiwan
Tel.:(886) 2-2624-2233 Fax.:(886) 2-2624-2266
www.fortunemed.com Mail:info@fortunemed.com