

SILICONE CLAMP JAW COVER

Provide cushions material of the clamp for hemostats, surgical clamps, and forceps.

SILICONE CLAMP JAW COVER


Fortune® Silicone Clamp Jaw Cover:

Silicone clamp jaw cover is used to provide the vessel or tissue with atraumatic clamping or occlusion of the vessel during surgery.

Suitable for surgeries using various clamps such as hemostatic forceps and surgical clamps.

Feature:


- It's made of medical-grade silicone, radiopaque, and available in 3 colors and 3 widths. The silicone elasticity can be used for most clamp sizes.
- The silicone clamp jaw cover can be trimmed to the desired length to fit the jaws of various-sized surgical instruments, and the serrated edge interface provides smooth, steady, non-slipping and firm grip.
- The shape of the lumen conforms with the shape of clamps or forceps, preventing the silicone clamp jaw cover from twisting or slipping.

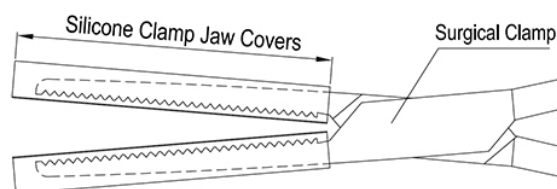

SPECIFICATIONS

REF. No.	Size	Length	Color
9002-SR250	2.5 mm	127 mm (5 inch)	Red
9002-MW410	4.1 mm	127 mm (5 inch)	White

REF. No.	Size	Length	Color
9002-LB570	5.7 mm	127 mm (5 inch)	Blue

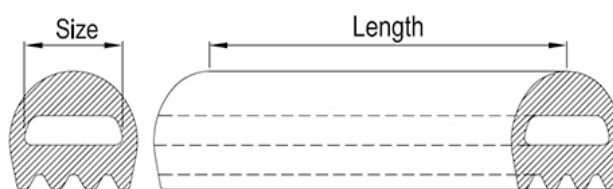

9002-XXXXX



6 F., No. 29, Sec. 2, Jhongjheng E. Rd.,
Danshuei Dist, New Taipei City. 251 Taiwan
Tel.:(886) 2-2624-2233 Fax.:(886) 2-2624-2266
www.fortunemed.com Mail:info@fortunemed.com